

medDispense® CT series

Expand your dispensing capacity



Overview

The medDispense® CT series automated dispensing cabinets allow you to expand your dispensing capacity while combining the security of single-item control and bulk supply storage.

- ✓ Accommodate both linear and matrix storage drawers combined with bulk storage
- ✓ Choose from five different cabinet sizes, three drawer options and secure bulk storage with a maximum capacity of 1,400 items
- ✓ Select the appropriate dispensing security for each item in the cabinet
- ✓ Powered by medLogic™ to deliver increased efficiency and flexibility for your team and a safer environment for your patients

Features

- Provides expanded capacity while combining the security of single item control and bulk supply storage
- Integrates with C and F series medDispense® automated dispensing cabinets
- Linear and matrix drawer options to choose from based on security level and item size
- Microban® antimicrobial infused surfaces
- Ability to easily combine bulk items with bins, and multiple drawer options in one cabinet
- The smallest unit fits on a countertop for areas with lower medication capacity needs
- Comprehensive 5- year warranty

Recommended for

- Psychiatric
- Critical access hospitals
- Simulation labs
- Long-term acute care

Other unit sizes in this series



Compatible with



C series



F series

medDispense® CT series

Expand your dispensing capacity

Accessories




External return bin



Overbridge



Bulk storage bins

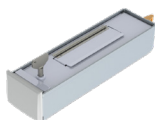

Return drawer

Specifications

Max unit capacity	112, 700, 800, or 1,400
Width	686 mm
Depth	559 mm (compact unit) or 711 mm (full depth unit)
Height	356 mm (112-capacity unit), 1,118 mm (700- and 800-capacity units) or 1,981 mm (1,400-capacity unit)
Weight	56-63 kgs (112-capacity unit), 143-206 kgs (700-capacity unit), 205-328 kgs (800-capacity unit), or 249-377 kgs (1,400-capacity unit)
Max number of drawers	12, 30, or 60 (dependent on unit capacity)
Material construction	Steel frame with polymer drawers and doors
Power	120VAC @ 60Hz

Software

- All medDispense® automated dispensing cabinets utilize the medLogic™ software platform
- The medLogic™ software platform delivers increased efficiency and flexibility for your team and a safer environment for your patients.
- **Efficiency:**
 - Single medication dispensing across multiple medDispense® automated dispensing cabinets.
 - Restock faster with barcode stocking.
 - FIFO stocking workflow ensures oldest stocked items are dispensed first, which helps reduced expired medication expenses.
- **Safety:**
 - Fully supports the UK NHS ECD register.
 - On-demand stock report generation by item/medication.
 - Help reduce medication errors with our configurable medication search.
 - Increase medication traceability with lot number stocking.
- **Flexibility:**
 - Fractional dose stocking & dispensing helps reduce medication waste.
 - High definition 16:9 and laptop display with the PharmacyCenter™ and NurSync™ modules provides increased technology flexibility for the pharmacy team.
- medLogic™ secures and tracks both your medications and supplies while seamlessly integrating with any hospital information system.
- medLogic™ modules improve workflow efficiency, safety and visibility so you can save time and focus on delivering exceptional patient care.
- **NurSync™ Remote Medication Queuing:**
 - Enables nurses to create medication queues for their patients to help reduce medication retrieval times, create a more efficient workflow and spend more time with patients.
- **medGuard™ Controlled Drug Inventory Management:**
 - Helps close the loop of controlled drug medication delivery and monitors for potential diversion with controlled restocking access and software safeguards.
 - Discrepancy tracking ensures the retrieval quantity recorded in medGuard™ matches the quantity restocked in the medDispense® automated dispensing cabinet.
- **PharmacyCenter™ Reporting and System Administration:**
 - Enables pharmacy staff to manage their medications across multiple medDispense® automated dispensing cabinets from one easy-to-use dashboard and reporting system.
 - Manage and maintain your facility's automated dispensing cabinets from any PC.
 - Visibility and management of your inventory across all automated dispensing cabinets.
 - Report on key metrics that matter to your facility, such as: inventory, top usage, transactions, controlled medications, and more.

